


## PTS Passer

*Suture Passer*

The Parcus PTS Passer helps to reliably pass suture through thick rotator cuff tissue. The side loading jaw is designed to be simple and easy to use, providing surgeons with intuitive operation and consistent control. The active window captures and retrieves suture for a simple one-handed pass. The PTS Passer fits down an 6.5mm cannula.

### Features & Benefits

- Side loading jaw for consistent control.
- Active window for a simple one-handed pass.
- Can pass either Parcus Suture Tape or suture.

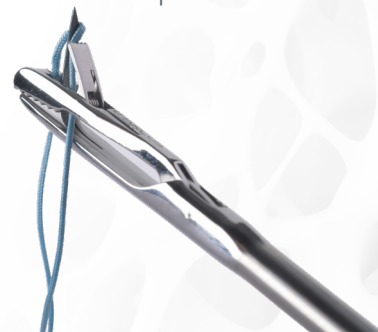

### Ordering Information

Part #	Description
11234	PTS Passer Needle, STERILE, Qty 1
20164	PTS Passer with Trap Door
20193	PTS Passer, Ratchet Handle with Trap Door

### Recommended Sutures and Tapes

Part #	Size	Description	Total Length (in)	Qty/Box	Working Loop
20032S	1.6mm	Parcus Braid, suture tape (wht/blu)	39	1	27
20034S	1.6mm	Parcus Braid, suture tape (wht/blk)	39	1	27
20042	2mm	Parcus Braid, suture tape (wht/blu)	39	6	7
20042S	2mm	Parcus Braid, suture tape (wht/blu)	39	1	7
20049	#2	Parcus Braid, suture (wht/blu & wht/blk)	36	12	-
20049S	#2	Parcus Braid, suture (wht/blu & wht/blk)	36	1	-

### Technique Tips

- Parcus Suture Tape is flat in the middle of the strand (tape) and round on either end (#2 suture).
- Use the rounded portion at the end of the strand of Parcus Suture Tape (the #2 suture) to pass through the tissue.
- Then slide the suture tape through so that the desired portion of the strand (the tape) is in contact with the tissue.



Parcus Medical has joined **Anika**

**PARCUS**

**Parcus Medical, LLC.**

6423 Parkland Dr, Sarasota, FL

1-941-755-7965 • [www.parcusmedical.com](http://www.parcusmedical.com)

[ParcusCS@anika.com](mailto:ParcusCS@anika.com)